


DIASONO-310

3/10 MHz **Multi-Frequency Modulated Ultrasonic Wave** and DIACELLPEEL's **Micro DermaBration feature** create a perfect combination of scaling, regeneration, sedation and lifting to provide safe and effective treatment and quick recuperative power.



Product Size / weight	Width 275 * Height 291 * Depth 350 / 9Kg
Power	100~240Vac, 50/60Hz, Power Consumption 800Wmax
Suction Pressure	MAX 750 mmHg
	Suction LEVEL : 10LEVEL
	Flow LEVEL : 10LEVEL
E/P	LEVEL : 10LEVEL
	LED : BLUE LED / When In Motion LED ON
Ultrasound	Frequency Selection 3MHz, 10MHz, 3/10MHz
	LEVEL : 10LEVEL
	Operation Time : MAX 30min
	MODE : 10 Program Mode / 2 User Mode
	LED : BLUE LED / When In Motion LED ON

SONO-310

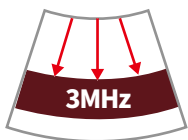


3MHz

It is transmitted up to 10mm depth of skin with 3 million vibrations per second, stimulating the dermis layers and delivering to the cell unit deep into the tissue for micro-massage effect.



10MHz

10 million vibrations per second, Nano Ultrasound focuses mainly 3mm depth of skin, the epidermis and dermis, to promote the synthesis of collagen and elastin, the major fibrous protein tissues that make up the skin.

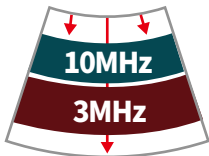

3/10MHz


3/10 MHz Multi-Frequency Modulated Ultrasound maximizes the benefits of Nano Ultrasound. As a technique, inhibits the proteolysis enzyme MMPs enzyme, increases the collagen-producing enzyme HSPs, and increases the GAGs to promote the growth and recovery of cells.

DIACELL PEEL



DIACELL PEEL

Using the exquisite whirlwind diamond tip made of natural diamond, the dead skin is carefully peeled like cutting jewel, and at the same time, tornado peeling is done with special 3 in 1 solution and powerful sound pressure function.

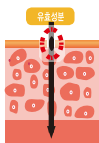

Electroporation

After replacing with an E/P tip, with electroporation function can make a thin hole in the skin and penetrate ampoule to desired area.



UV sterilization

It is equipped with a UV sterilization box that can safely disinfect accessories such as whirlwind diamond tips, so it has excellent hygiene function for skin safety.

Various adaptations and skin type programs are included

Can be used to specialize in acne-prone skin

Optimization of space and time utilization

Configurable for a variety of treatment

Easy procedure operation